

RP00158

Leader in Biomolecular Solutions for Life Science



Recombinant Human CD5 Protein

Catalog No.: RP00158 **Recombinant**

Sequence Information

Species	Gene ID	Swiss Prot
HEK293 cells	921	P06127

Tags

C-His

Synonyms

CD5; LEU1; T1; CD5 molecule; LEU1; T1

Product Information

Source	Purification
HEK293 cells	≥ 95 % as determined by SDS-PAGE.

Endotoxin

< 0.1 EU/μg of the protein by LAL method.

Formulation

Lyophilized from a 0.22 μm filtered solution of PBS, pH 7.4. Contact us for customized product form or formulation.

Reconstitution

Centrifuge the vial before opening. Reconstitute to a concentration of 0.1-0.5 mg/mL in sterile distilled water. Avoid vortex or vigorously pipetting the protein. For long term storage, it is recommended to add a carrier protein or stabilizer (e.g. 0.1% BSA, 5% HSA, 10% FBS or 5% Trehalose), and aliquot the reconstituted protein solution to minimize freeze-thaw cycles.

Contact



www.abclonal.com

Background

CD5, also known as Leu-1, Ly-1, and T1, is a 67 kDa transmembrane glycoprotein in the scavenger receptor superfamily. Members of this family are secreted or membrane-anchored proteins mainly found in cells associated with the immune system. This protein is a type-I transmembrane glycoprotein found on the surface of thymocytes, T lymphocytes and a subset of B lymphocytes. The encoded protein contains three SRCR domains and may act as a receptor to regulate T-cell proliferation.

Basic Information

Description

Recombinant Human CD5 Protein is produced by HEK293 cells expression system. The target protein is expressed with sequence (Arg25-Asn371) of human CD5 (Accession #NP_055022.2) fused with a 6×His tag at the C-terminus.

Bio-Activity

Shipping

The product is shipped at ambient temperature. Upon receipt, store it immediately at the temperature recommended below.

Operational Notes


For your safety and health, please wear a lab coat and disposable gloves for handling.

Storage

Store at -20°C. Store the lyophilized protein at -20°C to -80 °C up to 1 year from the date of receipt. After reconstitution, the protein solution is stable at -20°C for 3 months, at 2-8°C for up to 1 week. Avoid repeated freeze/thaw cycles.

* For your safety and health, please wear a lab coat and disposable gloves when handling.

Validation Data

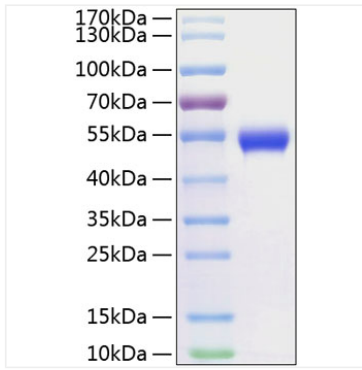

Recombinant Human CD5 Protein was determined by SDS-PAGE under reducing conditions with Coomassie Blue.