

A23772

Leader in Biomolecular Solutions for Life Science



# EV71 3C Rabbit mAb

Catalog No.: A23772

Recombinant

1 Publications

## Basic Information

---

### Observed MW

40kDa

### Calculated MW

20KDa

### Category

SMab Recombinant Monoclonal  
Antibody

### Applications

WB,ELISA

### Cross-Reactivity

Enterovirus

### CloneNo number

ARC59088

## Recommended Dilutions

---

**WB** 1:10000 - 1:170000

**ELISA** Recommended starting concentration is 1 µg/mL. Please optimize the concentration based on your specific assay requirements.

## Contact

---

 [www.abclonal.com](http://www.abclonal.com)

## Background

---

## Immunogen Information

---

### Gene ID

### Swiss Prot

### Immunogen

Recombinant protein (or fragment). This information is considered to be commercially sensitive.

### Synonyms

## Product Information

---

### Source

Rabbit

### Isotype

IgG

### Purification

Affinity purification


### Storage

Store at -20°C. Avoid freeze / thaw cycles.

Buffer: PBS containing 50% glycerol and 0.05% BSA, preserved with proclin300 or sodium azide (as specified on the Certificate of Analysis), pH 7.3.

## Validation Data

---

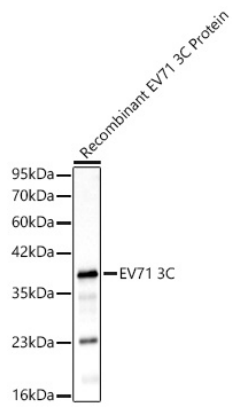

Western blot analysis lysates of recombinant EV71 3C Protein, using EV71 3C Rabbit mAb (A23772) at 1:160000 dilution.  
Secondary antibody: HRP-conjugated Goat anti-Rabbit IgG (H+L) (AS014) at 1:10000 dilution.  
Lysates/proteins: 5ng per lane.  
Blocking buffer: 3% nonfat dry milk in TBST.  
Detection: ECL Basic Kit (RM00020).  
Exposure time: 20s.