

A23087

Leader in Biomolecular Solutions for Life Science



Monkeypox virus (MPXV) E8L Rabbit mAb

Catalog No.: A23087

Recombinant

1 Publications

Basic Information

Observed MW

40-50kDa

Calculated MW

35kDa

Category

SMab Recombinant Monoclonal
Antibody

Applications

WB,ELISA

Cross-Reactivity

Monkeypox virus

CloneNo number

ARC59811

Recommended Dilutions

WB 1:2000 - 1:6000

ELISA Recommended starting concentration is 1 µg/mL. Please optimize the concentration based on your specific assay requirements.

Contact



www.abclonal.com

Background

Immunogen Information

Gene ID

929017

Swiss Prot

Q8V4Y0

Immunogen

Recombinant protein (or fragment). This information is considered to be commercially sensitive.

Synonyms

Product Information

Source

Rabbit

Isotype

IgG

Purification


Affinity purification

Storage

Store at -20°C. Avoid freeze / thaw cycles.

Buffer: PBS containing 50% glycerol and 0.05% BSA, preserved with proclin300 or sodium azide (as specified on the Certificate of Analysis), pH 7.3.

Validation Data

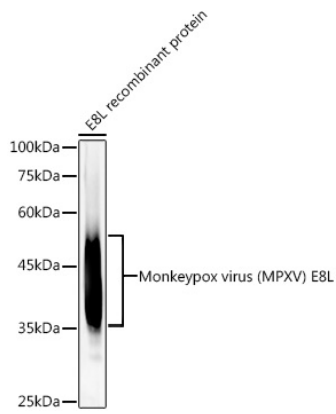

Western blot analysis of recombinant Monkeypox virus (MPXV) E8L protein, using Monkeypox virus (MPXV) E8L Rabbit mAb (A23087) at 1:5000 dilution. Secondary antibody: HRP-conjugated Goat anti-Rabbit IgG (H+L) (AS014) at 1:10000 dilution. Lysates/proteins: 10ng per lane. Blocking buffer: 3% nonfat dry milk in TBST. Detection: ECL Basic Kit (RM00020). Exposure time: 30s.