

AC028

Leader in Biomolecular Solutions for Life Science



HRP-conjugated β -Actin Rabbit mAb

Catalog No.: AC028

Recombinant

85 Publications

Basic Information

Observed MW

42 kDa

Calculated MW

42 kDa

Category

Monoclonal Antibody

Applications

WB,ELISA

Cross-Reactivity

Human,Mouse,Rat,Chicken,Zebrafish,Pig

CloneNo number

ARC5115-01

Conjugate

HRP

Recommended Dilutions

WB 1:50000-1:200000

ELISA Recommended starting concentration is 1 μ g/mL. Please optimize the concentration based on your specific assay requirements.

Contact



www.abclonal.com

Background

This gene encodes one of six different actin proteins. Actins are highly conserved proteins that are involved in cell motility, structure, integrity, and intercellular signaling. The encoded protein is a major constituent of the contractile apparatus and one of the two nonmuscle cytoskeletal actins that are ubiquitously expressed. Mutations in this gene cause Baraitser-Winter syndrome 1, which is characterized by intellectual disability with a distinctive facial appearance in human patients. Numerous pseudogenes of this gene have been identified throughout the human genome.

Immunogen Information

Gene ID

60

Swiss Prot

P60709

Immunogen

Recombinant protein (or fragment). This information is considered to be commercially sensitive.

Synonyms

BRWS1; PS1TP5BP1; HRP-conjugated β -Actin

Product Information

Source

Rabbit

Isotype

IgG

Purification


Affinity purification

Storage

Store at -20°C . Avoid freeze / thaw cycles.

Buffer: PBS containing 50% glycerol and 0.05% BSA, preserved with proclin300 or sodium azide (as specified on the Certificate of Analysis), pH 7.3.

Validation Data

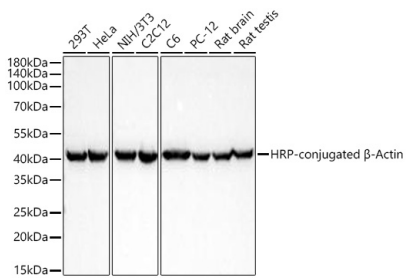

Western blot analysis of various lysates using HRP-conjugated β -Actin Rabbit mAb (AC028) at 1:200000 dilution incubated overnight at 4°C.
Secondary antibody: () at 1:10000 dilution.
Lysates/proteins: 25 μ g per lane.
Blocking buffer: 3% nonfat dry milk in TBST.
Detection: ECL Basic Kit (RM00020).
Exposure time: 10 s.