

AS076

Leader in Biomolecular Solutions for Life Science



www.abclonal.com

ABflo® 488-conjugated Goat anti-Mouse IgG (H+L)

Catalog No.: AS076

3 Publications

Basic Information

Observed MW

Calculated MW

Category

Secondary Antibody

Applications

IF/ICC,FC

Cross-Reactivity

Mouse

Conjugate

ABflo® 488. Ex:491nm. Em:516nm.

Background

Secondary antibodies are affinity-purified antibodies which will work with target-specific primary antibody in the detection, sorting or purification of its specified target. Secondary antibodies offer increased versatility enabling users to use many detection systems (e.g. HRP, AP, fluorescence). They can also provide greater sensitivity through signal amplification as multiple secondary antibodies. Most commonly, secondary antibodies are generated by immunizing the host animal (different from host species of primary antibody) with a pooled population of normal immunoglobulins from the host species of primary antibody and can be further purified and modified (i.e. antibody fragmentation, label conjugation, etc.) to ensure well-characterized specificity to corresponding normal immunoglobulins.

Recommended Dilutions

IF/ICC 1:50 - 1:200

FC 1:100 - 1:800

Immunogen Information

Gene ID

Swiss Prot

Immunogen

This information is considered to be commercially sensitive.

Synonyms

Contact



www.abclonal.com

Product Information

Source

Goat

Isotype

IgG

Purification


Affinity purification

Storage


Store at -20°C. Avoid freeze / thaw cycles.

Buffer: PBS with 0.025% Sodium Azide,0.75% BSA,50% glycerol,pH7.3.

Validation Data

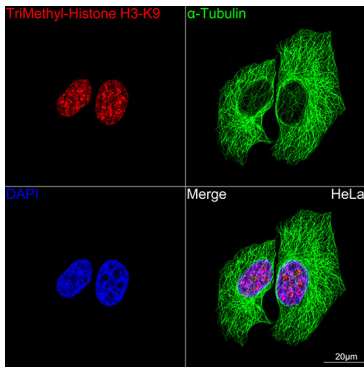

Confocal imaging of HeLa cells using TriMethyl-Histone H3-K9 Rabbit mAb (A22295,dilution 1:100)(Red). The cells were counterstained with α -Tubulin Mouse mAb (AC012,dilution 1:400) (Green). DAPI was used for nuclear staining (blue). Objective: 100x.Secondary antibody: Cy3 Goat Anti-Rabbit IgG (H+L) (AS007,dilution 1:100)(Red),ABflo® 488-conjugated Goat Anti-Mouse IgG (H+L)(AS076,dilution 1:200) (Green)

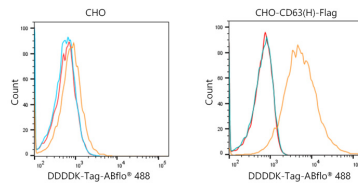

Flow cytometric analysis of Positive antibody (AE005) (2.5 μ g/mL) in various cells (orange) compare to Rabbit isotype control (blue) and non-staining control (Red).The secondary antibody used was ABflo® 488-conjugated AffiniPure Goat Anti-Mouse IgG (H+L) (AS076) at 1:100.