

AC038

Leader in Biomolecular Solutions for Life Science



β-Actin Rabbit mAb

Catalog No.: AC038

Recombinant

222 Publications

Basic Information

Observed MW

42 kDa

Calculated MW

42 kDa

Category

Monoclonal Antibody

Applications

WB,IHC-P,ELISA

Cross-Reactivity

Human,Mouse,Rat,Chicken,Zebrafish,Pig

CloneNo number

ARC5115-01

Recommended Dilutions

WB 1:10000 - 1:100000

IHC-P 1:400 - 1:5000

ELISA Recommended starting concentration is 1 µg/mL. Please optimize the concentration based on your specific assay requirements.

Contact

 www.abclonal.com

Background

This gene encodes one of six different actin proteins. Actins are highly conserved proteins that are involved in cell motility, structure, integrity, and intercellular signaling. The encoded protein is a major constituent of the contractile apparatus and one of the two nonmuscle cytoskeletal actins that are ubiquitously expressed. Mutations in this gene cause Baraitser-Winter syndrome 1, which is characterized by intellectual disability with a distinctive facial appearance in human patients. Numerous pseudogenes of this gene have been identified throughout the human genome.

Immunogen Information

Gene ID

60

Swiss Prot

P60709

Immunogen

Recombinant protein (or fragment). This information is considered to be commercially sensitive.

Synonyms

BRWS1; PS1TP5BP1; β-Actin

Product Information

Source

Rabbit

Isotype

IgG

Purification


Affinity purification

Storage


Store at -20°C. Avoid freeze / thaw cycles.

Buffer: PBS containing 50% glycerol and 0.05% BSA, preserved with proclin300 or sodium azide (as specified on the Certificate of Analysis), pH 7.3.


Validation Data

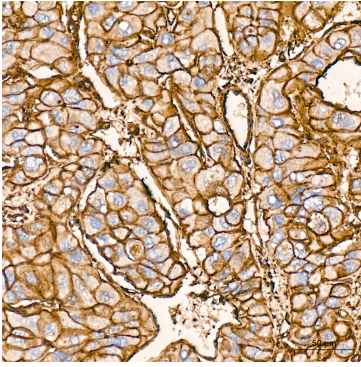


Immunohistochemistry analysis of paraffin-embedded Human liver cancer tissue using β -Actin Rabbit mAb (AC038) at a dilution of 1:400 (40x lens). High pressure antigen retrieval performed with 0.01M Citrate Buffer (pH 6.0) prior to IHC staining.

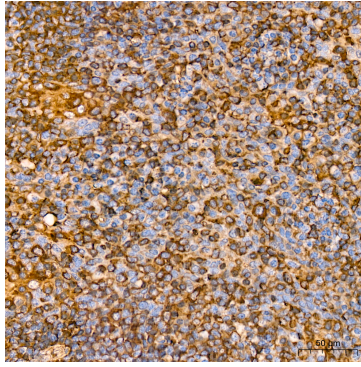

Immunohistochemistry analysis of paraffin-embedded Human tonsil tissue using β -Actin Rabbit mAb (AC038) at a dilution of 1:400 (40x lens). High pressure antigen retrieval performed with 0.01M Citrate Buffer (pH 6.0) prior to IHC staining.

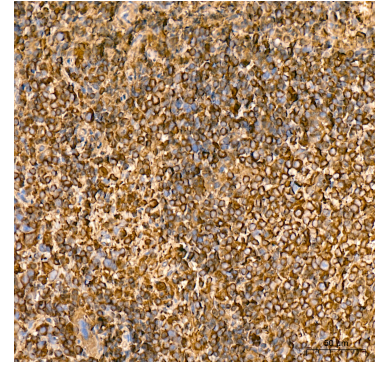

Immunohistochemistry analysis of paraffin-embedded Mouse spleen tissue using β -Actin Rabbit mAb (AC038) at a dilution of 1:400 (40x lens). High pressure antigen retrieval performed with 0.01M Citrate Buffer (pH 6.0) prior to IHC staining.

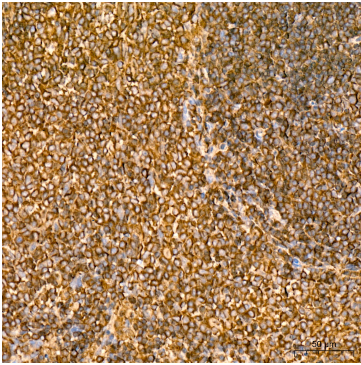

Immunohistochemistry analysis of paraffin-embedded Rat spleen_ tissue using β -Actin Rabbit mAb (AC038) at a dilution of 1:400 (40x lens). High pressure antigen retrieval performed with 0.01M Citrate Buffer (pH 6.0) prior to IHC staining.