

A1325

Leader in Biomolecular Solutions for Life Science



CXCL12 Rabbit pAb

Catalog No.: A1325 **9 Publications**

Basic Information

Observed MW

11kDa

Calculated MW

11kDa

Category

Polyclonal Antibody

Applications

WB,IHC-P,IF/ICC,ELISA

Cross-Reactivity

Human,Mouse,Rat

Background

This antimicrobial gene encodes a stromal cell-derived alpha chemokine member of the intercrine family. The encoded protein functions as the ligand for the G-protein coupled receptor, chemokine (C-X-C motif) receptor 4, and plays a role in many diverse cellular functions, including embryogenesis, immune surveillance, inflammation response, tissue homeostasis, and tumor growth and metastasis. Mutations in this gene are associated with resistance to human immunodeficiency virus type 1 infections. Multiple transcript variants encoding different isoforms have been found for this gene.

Recommended Dilutions

WB 1:500 - 1:1000

IHC-P 1:100 - 1:200

IF/ICC 1:100 - 1:200

ELISA Recommended starting concentration is 1 µg/mL. Please optimize the concentration based on your specific assay requirements.

Immunogen Information

Gene ID

6387

Swiss Prot

P48061

Immunogen

Synthetic peptide. This information is considered to be commercially sensitive.

Synonyms

IRH; PBSF; SDF1; TLSF; TPAR1; SCYB12; CXCL12

Product Information

Source

Rabbit

Isotype

IgG

Purification

Affinity purification

Storage

Store at -20°C. Avoid freeze / thaw cycles.


Buffer: PBS with 0.09% sodium azide,50% glycerol,pH7.3.

Contact



www.abclonal.com

Validation Data

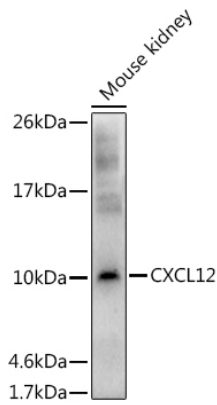

Western blot analysis of lysates from Mouse kidney, using CXCL12 Rabbit pAb (A1325) at 1:1000 dilution.


Secondary antibody: HRP-conjugated Goat anti-Rabbit IgG (H+L) (AS014) at 1:10000 dilution.

Lysates/proteins: 25µg per lane.


Blocking buffer: 3% nonfat dry milk in TBST.

Detection: ECL Enhanced Kit (RM00021).


Exposure time: 180s.

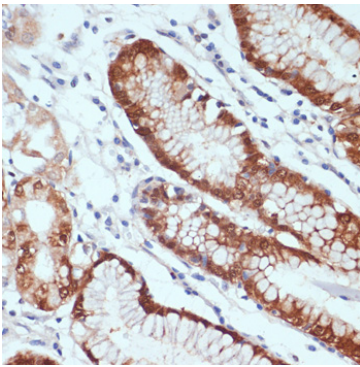


Immunohistochemistry analysis of paraffin-embedded Human stomach using CXCL12 Rabbit pAb (A1325) at dilution of 1:100 (40x lens). Microwave antigen retrieval performed with 0.01M PBS Buffer (pH 7.2) prior to IHC staining.

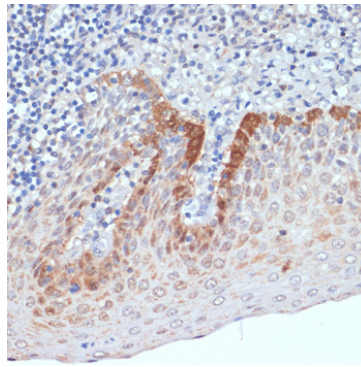

Immunohistochemistry analysis of paraffin-embedded Human tonsil using CXCL12 Rabbit pAb (A1325) at dilution of 1:100 (40x lens). Microwave antigen retrieval performed with 0.01M PBS Buffer (pH 7.2) prior to IHC staining.

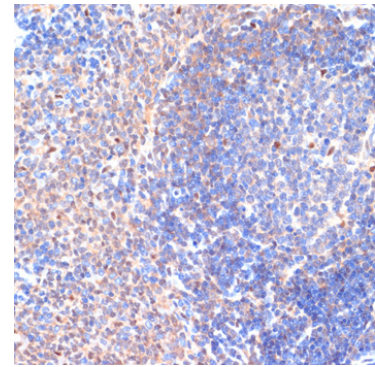

Immunohistochemistry analysis of paraffin-embedded Rat spleen using CXCL12 Rabbit pAb (A1325) at dilution of 1:100 (40x lens). Microwave antigen retrieval performed with 0.01M PBS Buffer (pH 7.2) prior to IHC staining.

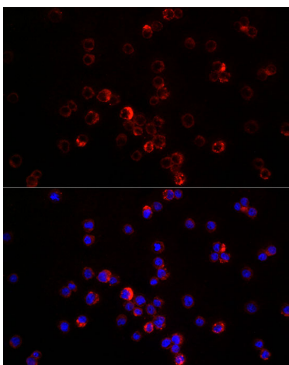

Immunofluorescence analysis of THP-1 cells using CXCL12 Rabbit pAb (A1325) at dilution of 1:200. Secondary antibody: Cy3-conjugated Goat anti-Rabbit IgG (H+L) (AS007) at 1:500 dilution. Blue: DAPI for nuclear staining.