

A18015

Leader in Biomolecular Solutions for Life Science



[KO Validated] Hsp90 β Rabbit pAb

Catalog No.: A18015

KO Validated

Basic Information

Observed MW

80kDa

Calculated MW

83kDa

Category

Polyclonal Antibody

Applications

WB,IF/ICC,ELISA

Cross-Reactivity

Human,Mouse,Rat,Monkey

Background

This gene encodes a member of the heat shock protein 90 family; these proteins are involved in signal transduction, protein folding and degradation and morphological evolution. This gene encodes the constitutive form of the cytosolic 90 kDa heat-shock protein and is thought to play a role in gastric apoptosis and inflammation. Alternative splicing results in multiple transcript variants. Pseudogenes have been identified on multiple chromosomes.

Recommended Dilutions

WB 1:500 - 1:2000

IF/ICC 1:50 - 1:200

ELISA Recommended starting concentration is 1 μ g/mL. Please optimize the concentration based on your specific assay requirements.

Immunogen Information

Gene ID

3326

Swiss Prot

P08238

Immunogen

Recombinant protein (or fragment). This information is considered to be commercially sensitive.

Synonyms

HSP84; HSPC2; HSPCB; D6S182; HSP90B; O β

Contact

 www.abclonal.com

Product Information

Source

Rabbit

Isotype

IgG

Purification


Affinity purification

Storage

Store at -20°C. Avoid freeze / thaw cycles.

Buffer: PBS with 0.01% thimerosal, 50% glycerol, pH7.3.

Validation Data

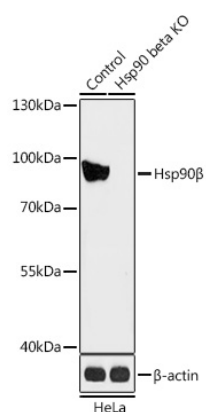

Western blot analysis of lysates from wild type (WT) and HSP90AB1 knockout (KO) HeLa cells, using [KO Validated] Hsp90β Rabbit pAb (A18015) at 1:1000 dilution.

Secondary antibody: HRP-conjugated Goat anti-Rabbit IgG (H+L) (AS014) at 1:10000 dilution.

Lysates/proteins: 25μg per lane.

Blocking buffer: 3% nonfat dry milk in TBST.

Detection: ECL Basic Kit (RM00020).

Exposure time: 90s.