


# SensiTek Anti-Polyvalent

**Description:** Polyclonal secondary antibody conjugated to biotin for 3-step immunohistochemistry protocols. Formulated to provide optimal staining with an incubation for 15-20 minutes. May be used with automated systems, reagent jars, and manual dropping/pipetting.

Species of Origin: Goat  
 Antigen Specificity: Anti-Rabbit, Mouse, Rat, Rabbit, Guinea Pig IgG (H+L)  
 Preadsorbed Against: Human

**Uses/Limitations:** Not to be taken internally.  
 For In-Vitro Diagnostic use only.  
 Histological applications.  
 Do not use if reagent becomes cloudy.  
 Do not use past expiration date.  
 Use caution when handling reagent.  
 Non-Sterile.

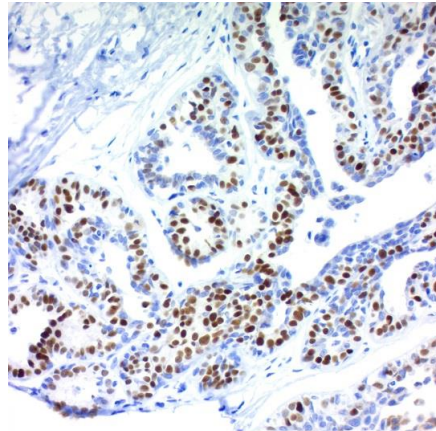

Human Breast Carcinoma stained with DAB and SensiTek Anti-Polyvalent within an IHC protocol.

**Control Tissue:** Any well-fixed tissue.

**Availability/Contents:**

<u>Item #</u>	<u>Volume</u>
ABF008	8 ml
ABF015	15 ml
ABF125	125 ml
ABF500	500 ml
ABF999	1000 ml

**Storage:** Store at 2-8°C. Product is stable for 18 months from date of manufacture.

**Procedure:**

*Allow reagents to come to room temperature before use.*

1. Deparaffinize and rehydrate tissue section.
2. If needed, incubate slide in hydrogen peroxide for 10-15 minutes to reduce nonspecific background staining due to endogenous peroxidase.
3. Wash 2 times in buffer.
4. If required, incubate tissue in digestive enzyme.
5. Wash 4 times in buffer.
6. Place slide in protein block (ScyTek's "Superblock") and incubate 5-10 minutes at room temperature to block nonspecific background staining. Note: Do not exceed 10 minutes or there may be a reduction in desired stain.

Storage: 2°C  8°C

P.O. Box 3286 - Logan, Utah 84323, U.S.A. - Tel. (800) 729-8350 - Fax (435) 755-0015 - [www.scytek.com](http://www.scytek.com)

7. Wash 1 time in buffer.
8. Apply primary antibody and incubate according to manufacturer's protocol.
9. Wash 4 times in buffer.
10. Place slide in SensiTek Anti-Polyvalent, and incubate for 15-20 minutes at room temperature.
11. Wash 4 times in buffer.
12. Place slide in enzyme label, and incubate per instructions
13. Rinse 4 times in buffer.
14. Place slide in appropriate chromogenic substrate and incubate until desired reaction is achieved.
15. Counterstain and coverslip.

### Troubleshooting Guide

#### **Overstaining:**

1. Concentration of the primary antibody was too high or the incubation time was too long.
2. Temperature during incubation was too high.
3. Incubation time with link antibody or streptavidin/enzyme label was too long.

#### **Nonspecific Background Staining:**

1. Rinsing between steps was inadequate.
2. Tissue was allowed to dry with reagents on.
3. Folds in tissue trapped reagents.
4. Tissue contains endogenous enzyme.
5. Tissue contains endogenous biotin.
6. Antigen migrated in tissue.
7. Excessive tissue adhesive on slides.
8. Inadequate blocking with protein block.

#### **Weak Staining:**


1. Primary antibody concentration was too low or incubation time was too short.
2. Reagents are past their expiration date.
3. Reagent is reaching the end of its useful life.
4. Counterstain or mounting media were incompatible and dissolved the chromogen reaction product.
5. Room temperature was excessively cool.
6. The primary antibody does not recognize an antigen that survives fixation and embedding in high enough amounts.
7. Excessive incubation with protein block (Super Block or normal serum).

#### **No Staining:**

1. Steps were inadvertently left out.
2. There is no relevant antigen in the tissue.
3. Chromogenic substrate does not match enzyme label.
4. One or more components have been inactivated.

#### **References:**

Storage: 2°C  8°C



ScyTek Laboratories, Inc.  
205 South 600 West  
Logan, UT 84321  
U.S.A.



EC REP  
Emergo Europe  
Westervoortsedijk 60  
6827 AT Arnhem, The Netherlands

Instructions For Use  
ABF-IFU


Rev. Date: July 2, 2019

Revision: 4

Page 3 of 3

P.O. Box 3286 - Logan, Utah 84323, U.S.A. - Tel. (800) 729-8350 - Fax (435) 755-0015 - [www.scytek.com](http://www.scytek.com)

1. Günes-Bayir, A., Kocyigit, A., Güler, E. M., Bilgin, M. G., Ergün, İ. S., & Dadak, A. (2018). Effects of carvedilol on human fibroblast (WS-1) and gastric adenocarcinoma (AGS) cells in vitro and on Wistar rats in vivo. *Molecular and Cellular Biochemistry*, 448(1), 237–249. <https://doi.org/10.1007/s11010-018-3329-5>
2. Koç, M., Kumral, Z. N. Ö., Özkan, N., Memi, G., Kaçar, Ö., Bilsel, S., ... Yeğen, B. Ç. (2014). Obestatin improves ischemia/reperfusion-induced renal injury in rats via its antioxidant and anti-apoptotic effects: Role of the nitric oxide. *Peptides*, 60, 23–31. <https://doi.org/10.1016/j.peptides.2014.07.019>
3. Uzar, E., Acar, A., Evliyaoğlu, O., Fırat, U., Kamasak, K., Göçmez, C., ... İlhan, A. (2012). The anti-oxidant and anti-apoptotic effects of nebivolol and zofenopril in a model of cerebral ischemia/reperfusion in rats. *Progress in Neuro-Psychopharmacology and Biological Psychiatry*, 36(1), 22–28. <https://doi.org/10.1016/j.pnpbp.2011.08.011>
4. Yasoda, A., Ogawa, Y., Suda, M., Tamura, N., Mori, K., Sakuma, Y., ... Nakao, K. (1998). Natriuretic Peptide Regulation of Endochondral Ossification EVIDENCE FOR POSSIBLE ROLES OF THE C-TYPE NATRIURETIC PEPTIDE/GUANYLYL CYCLASE-B PATHWAY. *Journal of Biological Chemistry*, 273(19), 11695–11700. <https://doi.org/10.1074/jbc.273.19.11695>
5. Yoshida, T., Matsumoto, E.-I., Hanamura, N., Kalembeiyi, I., Katsuta, K., Ishihara, A., & Sakakura, T. (1997). Co-expression of tenascin and fibronectin in epithelial and stromal cells of benign lesions and ductal carcinomas in the human breast. *The Journal of Pathology*, 182(4), 421–428. [https://doi.org/10.1002/\(SICI\)1096-9896\(199708\)182:4<421::AID-PATH886>3.0.CO;2-U](https://doi.org/10.1002/(SICI)1096-9896(199708)182:4<421::AID-PATH886>3.0.CO;2-U)

Storage:  
2°C  8°C ScyTek Laboratories, Inc.  
205 South 600 West  
Logan, UT 84321  
U.S.A.CE  Emergo Europe  
Westervoortsedijk 60  
6827 AT Arnhem, The Netherlands