

Distributed by:

CliniSciences Group

## CONTENIS

01

16

## Nucleic Acid Purification

MagicPure® Blood Genomic DNA Kit
MagicPure® Simple Viral DNA/RNA Kit
MagicPure® Cell-Free DNA Kit II
MagicPure® Stool and Soil Genomic DNA Kit
MagicPure® Soil Genomic DNA Kit
MagicPure® FFPE Tissue Genomic DNA Kit
MagicPure® Fungi Genomic DNA Kit
EasyPure® Stool Genomic DNA Kit
EasyPure® Stool Genomic DNA Kit
EasyPure® Buccal Swab Genomic DNA Kit
EasyPure® Viral DNA/RNA Kit
MagicPure® Host Cell Residual DNA Kit

09

## Common Used Enzymes

PerfectStart® Taq DNA Polymerase
TransStart® FastPfu DNA Polymerase
TransStart® FastPfu Fly DNA Polymerase
Bst II DNA Polymerase
Bst III DNA Polymerase
T4 DNA Ligase for NGS
DNA Polymerase I Klenow Fragment
T4 DNA Polymerase
T4 Polymerase
T4 Polynucleotide Kinase
Universal Nuclease (GMP Grade)
Uracil-DNA Glycosylase (Low Temperature)

## RT-PCR

TransScript® IV One-Step gDNA Removal and cDNA Synthesis SuperMix

TransScript® One-Step gDNA Removal and cDNA Synthesis SuperMix

TransScript®miRNA First-Strand cDNA Synthesis SuperMix

19

## qPCR and qRT-PCR

TransScript® II Multiplex Probe One-Step qRT-PCR SuperMix UDG

PerfectStart® II Probe qPCR SuperMix UDG PerfectStart® Fast Green qPCR SuperMix

**22** 

## In Vitro mRNA Synthesis

T7 High Efficiency Transcription Kit mRNA Capping Kit mRNA Poly(A) Tailing Kit


# TRANSGEN Nucleic Acid Purification



MagicPure® Blood Genomic DNA Kit (EC101)

- Simple and fast, no centrifugation required.
- High yield and purity.
- Suitable for fresh, frozen or anticoagulated blood.

## **Comparison with Competitive Products**

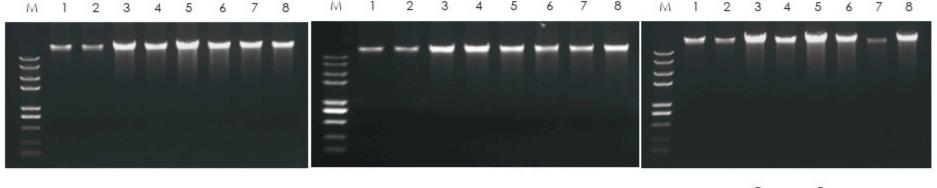

M: Trans2K<sup>®</sup> Plus II DNA Marker 1, 2: Purified nucleic acids from blood sample #A

3, 4: Purified nucleic acids from

blood sample #B 5, 6: Purified nucleic acids from

blood sample #C 7, 8: Purified nucleic acids from blood sample #D


TransGen Company A Company B


Genomic DNA was purified from 200 µl of human whole blood using purification kits from TransGen, Company A and Company B respectively. The purified nucleic acids were analyzed by 1.5% agarose gel electrophoresis.

## MagicPure® Simple Viral DNA/RNA Kit (EC311)

- Simple, no centrifugation required.
- Fast, high yield.
- Silica-based magnetic beads to specifically adsorb and purify viral DNA/RNA. No proteinase K treatment needed. It is suitable for high-throughput nucleic acid purification.
- Suitable for plasma, serum, whole blood, tissue homogenate, cell-free body fluid, nasopharyngeal or oropharyngeal aspirate/ wash, bronchoalveolar lavage fluid (BALF), tracheal aspirate, sputum.



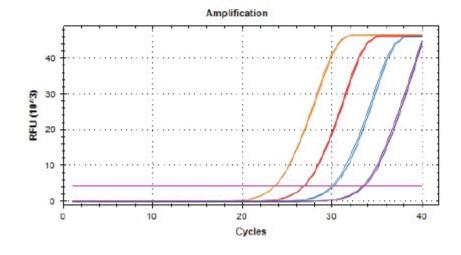

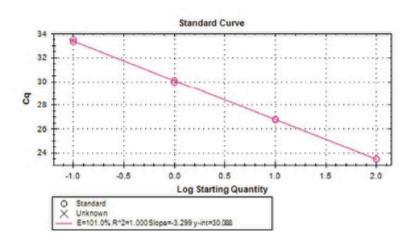


MagicPure Simple Viral DNA/RNA Kits were used to extract nucleic acids from 10<sup>1</sup> 10<sup>2</sup>, 10<sup>3</sup>, 10<sup>4</sup>-fold dilutions of porcine pseudorabies virus (PRV). The extracted nucleic acids were analyzed by qPCR.


## 

## MagicPure® Cell-Free DNA Kit II (EC211)


- Simple and fast, no centrifugation required.
- High yield and purity.
- Suitable for human serum or plasma.

Main cfDNA fragment concentration measured by Agilent 2100

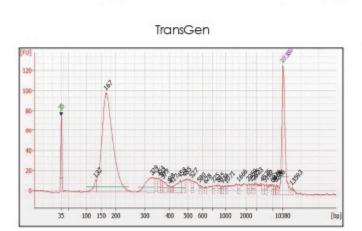


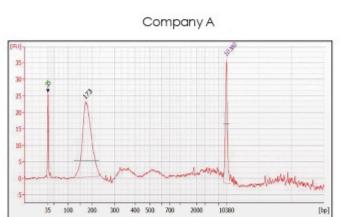


Library concentration measured by Qubit

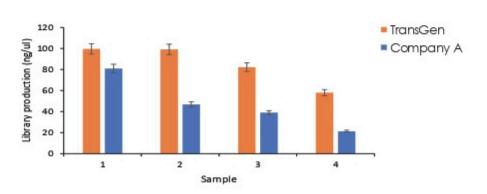

Library concentration measured by Agilent 4200

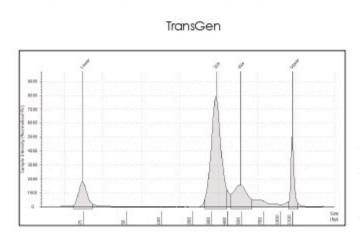

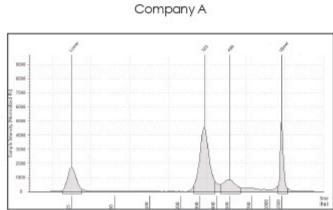


## MagicPure® Stool and Soil Genomic DNA Kit (EC801)

- · Fast extraction.
- High yield and high purity, combined optimized buffer with magnetic beads to achieve high yield and remove inhibitors.
- Suitable for fresh, frozen solid or liquid stool and soil.


## Comparison with Competitive Product

## Soil Samples

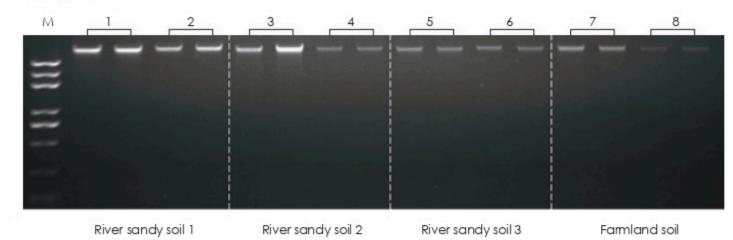


M: Trans2K<sup>®</sup> Plus DNA Marker 1, 3, 5, 7: TransGen 2, 4, 6, 8: Company A

Genomic DNA was extracted from 250 mg of different types of soil samples using reagents from TransGen and Company A respectively. The extracted nucleic acids were analyzed by agarose gel electrophoresis.

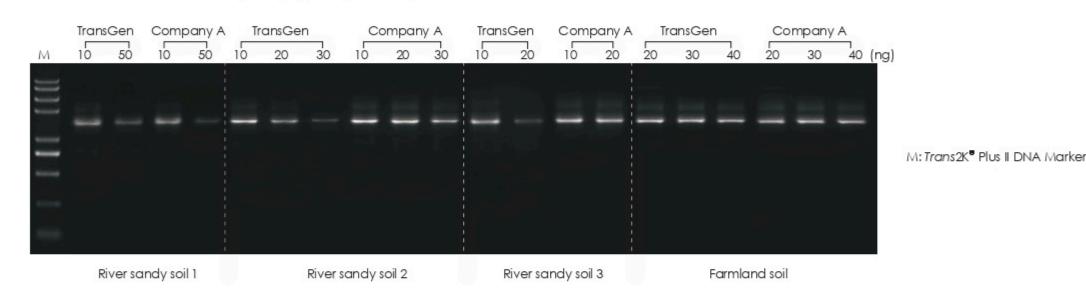

Genomic DNA was extracted from different types of soil samples using reagents from TransGen and Company A respectively. Extracted genomic DNA was used as template for 16s rRNA gene amplification by PCR. The amplified nucleic acids were analyzed by agarose gel electrophoresis.


## Stool Samples

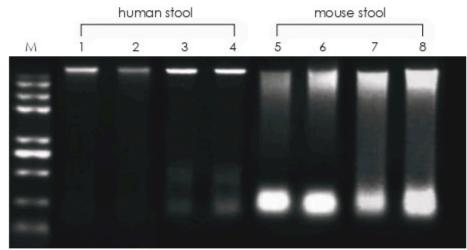

M: Trans2K<sup>®</sup> Plus II DNA Marker 1, 2, 5, 6: TransGen

3, 4, 7, 8: Company A

Genomic DNA was extracted from 200 mg of mouse and 200 µl of human stool using reagents from TransGen and Company A respectively. The extracted nucleic acids were analyzed by agarose gel electrophoresis.

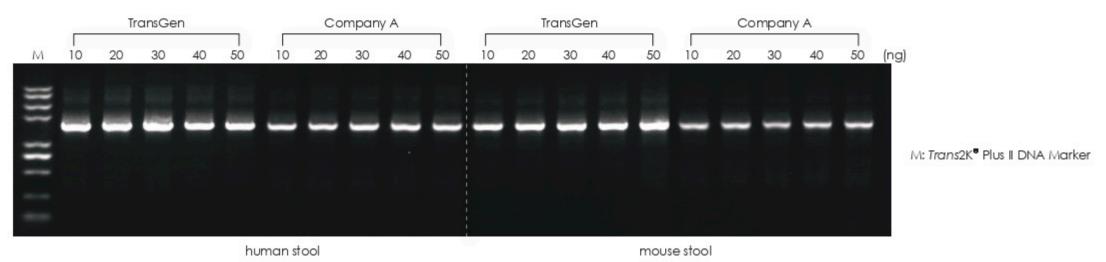

Genomic DNA was extracted from human and mouse stool samples using reagents from TransGen and Company A respectively. Extracted genomic DNA was used as template for 16s rRNA gene amplification by PCR. The amplified nucleic acids were analyzed by agarose gel electrophoresis.

MagicPure® Soil Genomic DNA Kit

(EC802)


- High purity, the magnetic bead method is used to specifically adsorb DNA, which can efficiently remove inhibitors in the sample.
- High extraction efficiency: suitable for efficient extraction of soil DNA with low microbial content.
- Wide range of application: adding Lysis Enhancer components to improve the lysis ability of special soil samples, suitable for the extraction of various types of soil samples.




TransGen

Company B

## Results

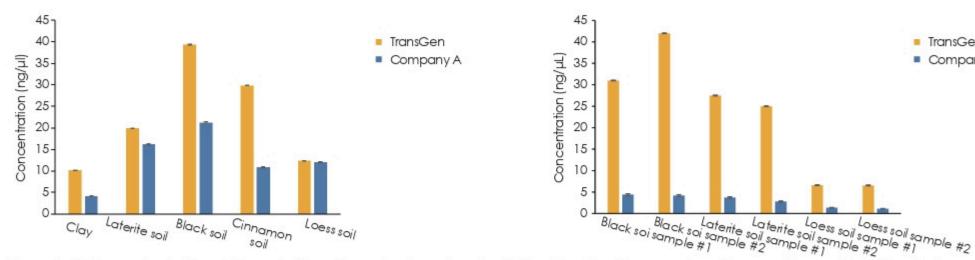


Genomic DNA was extracted from 250 mg of different types of soil samples using kits from TransGen, Company A and Company B respectively. Extracted genomic DNA was used as template for 16s rRNA gene amplification by PCR. The amplification nucleic acids were analyzed by agarose gel electrophoresis.

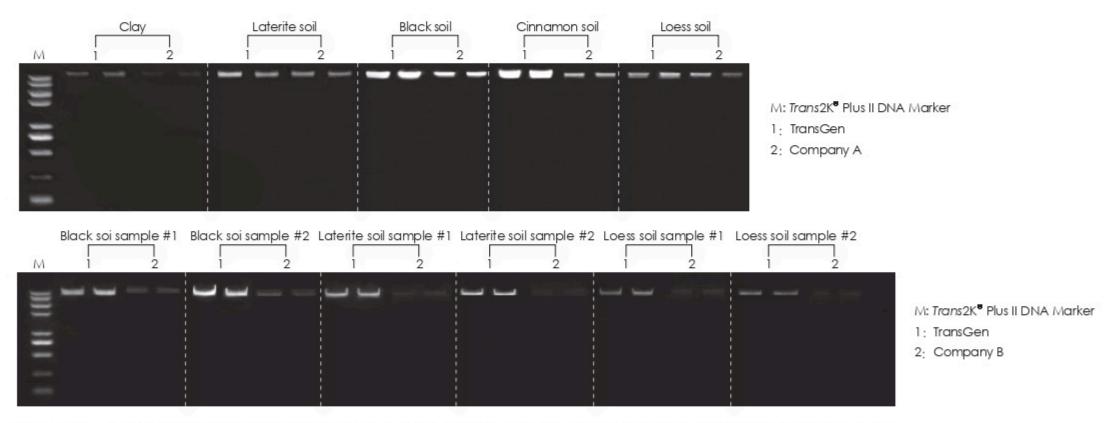

Genomic DNA was extracted from 250 mg of different types of soil samples using kits from TransGen, Company A and Company B respectively. The extracted gDNA was analyzed by electrophoresis.

## No Inhibition of Downstream PCR Amplification

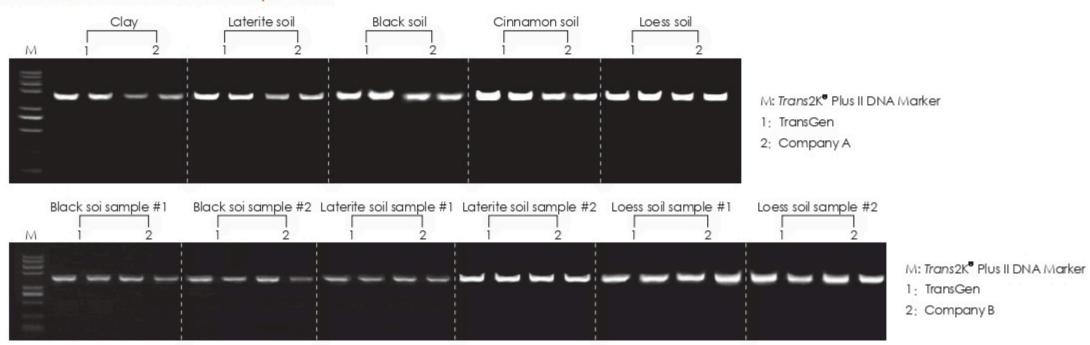

Different types of soil sample gDNA, extracted by reagents from TransGen, Company A and Company B, were used as templates to amplify 16S rDNA by PCR. The amplified nucleic acids were analyzed by agarose gel electrophoresis.


## MagicPure® FFPE Tissue Genomic DNA Kit (EC701)

- Fast, high yield.
- High purity enabled by buffer optimized for paraffin-embedded tissue samples and magnetic beads allowing for high-efficiency and specific DNA adsorption, and effective removal of protein, salt, and other inhibitors.
- Suitable for fresh formalin-fixed and paraffinembedded tissues and paraffin sections.

## MagicPure® Fungi Genomic DNA Kit (EC103)


- · Simple and fast.
- High yield and purity, with unique lysis buffer.
- Ideal for efficient extraction of genomic DNA from a variety of fungi such as Fusarium, yeast, anthrax and others.

## **Comparsion with Competitive Product**

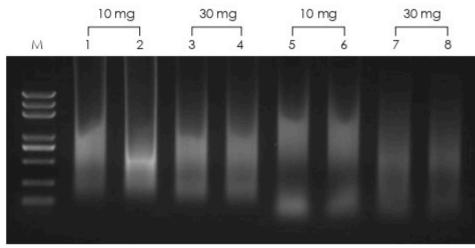

M: Trans2K<sup>®</sup> Plus DNA Marker 1, 2, 3, 4: TransGen 5, 6, 7, 8: Company A

Genomic DNA was extracted from different amounts of human paraffin samples using reagnets from TransGen and Company A respectively. The extracted nucleic acids were analyzed by agarose gel electrophoresis

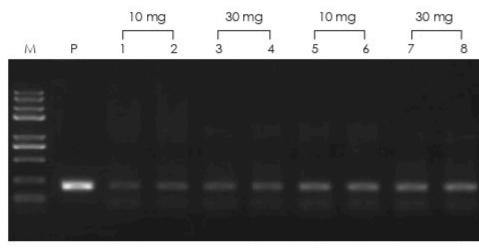


M: Trans2K<sup>®</sup> Plus II DNA Marker P: Positive Control (Human Blood gDNA)

1, 2, 3, 4: TransGen

5, 6, 7, 8: Company A


Genomic DNA was extracted from different amounts of human paraffin samples using reagents from TransGen and Company A respectively. The extracted gDNAs (20 ng) were used as templates for PCR. The amplified nucleic acids were analyzed by agarose gel electrophoresis.

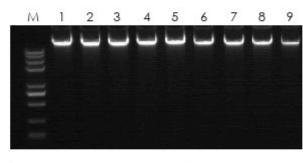

M: Trans2K<sup>®</sup> Plus II DNA Marker 1-3: Fusarium

4-6:Yeast 7-9:Anthrax

Genomic DNA was extracted from Fusarium, yeast, and anthrax using MagicPure® Fungi Genomic DNA Kit., The extracted nucleic acids were analyzed by agarose gel electrophoresis.

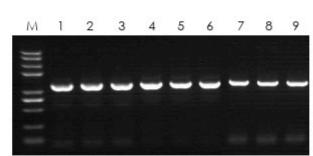

M: Trans2K<sup>®</sup> Plus II DNA Marker 1-3:Fusarium (NPS6)

4-6:Yeast(β-Actin)


7-9:Anthrax(RPOB)

Genomic DNA was extracted from Fusarium, yeast, and anthrax using MagicPure Fungi Genomic DNA Kit. The extracted gDNAs were used as templates for PCR. The amplified nucleic acids were analyzed by agarose gel electrophoresis.

EasyPure® Stool Genomic DNA Kit (EE301)


- Fast extraction, good DNA integrity, high yield.
- High purity, buffer optimized for stool samples, combined with magnetic beads which bind to DNA efficiently and specifically, can remove the inhibitors in the samples effectively.
- Suitable for fresh, frozen solid or liquid stool.

## **Comparison with Competitive Products**

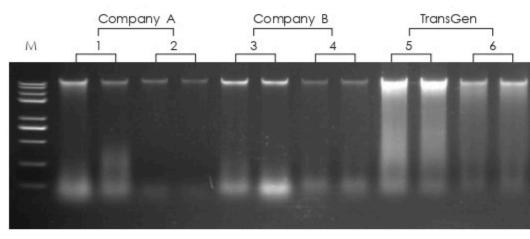

M: Trans2K<sup>®</sup> Plus II DNA Marker
1, 3, 5: Mouse stool
2, 4, 6: Human stool

Genomic DNA was extracted from stool samples of mouse and human using reagents from TransGen, Company A and Company B respectively. The extracted nucleic acids were analyzed by agarose gel electrophoresis.

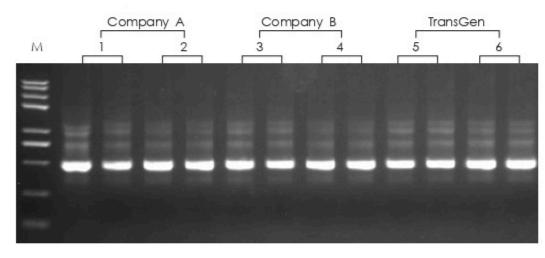

M: Trans2K® Plus II DNA Marker

1, 3, 5: Mouse stool


2, 4, 6: Human stool

Genomic DNA was extracted from mouse and human stool samples using reagents from TransGen, Company A and Company B respectively. Genomic DNA for 16s rRNA gene was amplified. The amplification products were analyzed by agarose gel electrophoresis.

## EasyPure® Buccal Swab Genomic DNA Kit (EE201)

- · Fast extraction, high yield.
- High purity, enabled by spin column with specific adsorption to DNA, effectively removing impurities as protein and salts, and others.

## Comparsion with Competitive Product

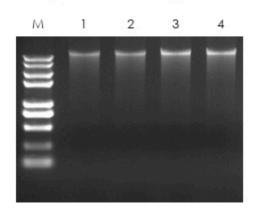

M: Trans2K® Plus II DNA Marker

1, 2: Company A

3, 4: TransGen

Genomic DNA was extracted from buccal swabs using reagents of TransGen and Company A respectively. The extracted nucleic acids were analyzed by 1.0% agarose gel electrophoresis.




## EasyPure® Viral DNA/RNA Kit (ER201)


- · Easy procedures and high yield
- Use a unique lysis buffer to lyse virus and silicabased spin column to specifically adsorb DNA.
- Suitable for serum, plasma, whole blood, tissue homogenates, cell-free body fluids, nasopharyngeal or oropharyngeal aspirates or washes, alveolar lavage, tracheal aspirates and sputum, nasopharyngeal or oropharyngeal swabs, culture animal cell supernatant.

## MagicPure® Host Cell Residual DNA Kit (EH201)

- Use a unique lysate to release a small amount of DNA, and use silica-based magnetic beads to specifically adsorb and purify residual DNA from host cells.
- Simple operation, no need for centrifugation and fast extraction.
- High extraction yield, for the extraction of trace DNA fragments from host cells (CHO, Vero, NSO, MDCK, etc.) from ≤200 µL biological samples.
- High purity, suitable for experiments such as PCR, RT-PCR, qPCR and qRT-PCR.
- It is suitable for high-throughput magnetic rod-type nucleic acid extractors.


## Comparsion with Competitive Product

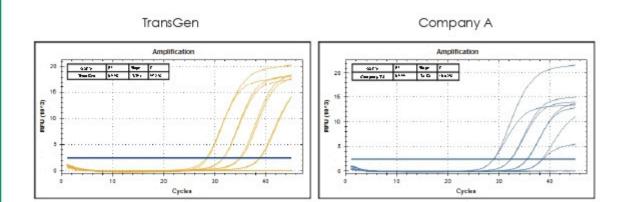

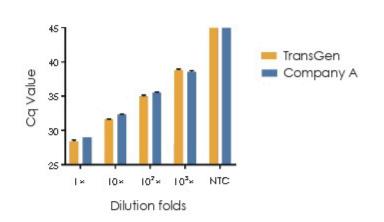

DNA was extracted from African Swine Fever Virus (ASFV) infected swine serum using reagents from TransGen and Company A. The extracted gDNAs were used as templates for qPCR.

## **Comparsion with Competitive Product**

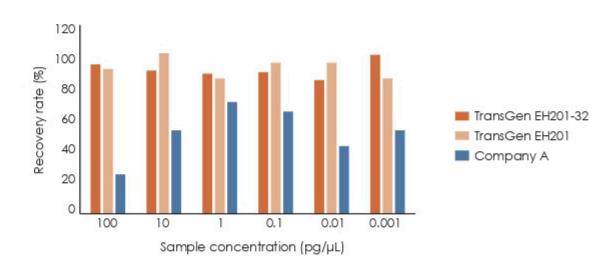

A 10-fold gradient dilution of Vero fragmented DNA (sonicated DNA) containing a high concentration of protein was used to extract DNA with MagicPure® Host Cell Residual DNA Kit and with reagent from Company A. The extracted DNA concentrations were measured. TransGen EH201-32 is used with a 32-channel automatic nucleic acid extractor.

## TRANSGEN

## Common Used Enzymes




## PerfectStart® Taq DNA Polymerase (AP401)

- · Hot start and high specificity.
- · High sensitivity.
- · High amplification efficiency.
- Low E. coli genomic DNA residue.
- Suitable for low copy number templates, complex templates, GC/AT-rich templates, real-time quantitative PCR (qPCR), and multiplex PCR and single nucleotide polymorphisms (SNPs)



## **Comparsion with Competitive Product**

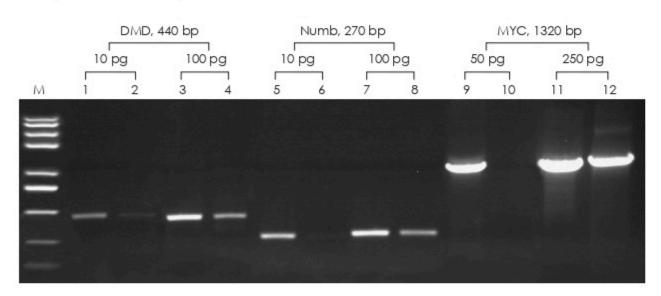


M: Trans2K® Plus II DNA Marker 1, 3, 5, 7, 9, 11: TransGen 2, 4, 6, 8, 10, 12: Company A

Genomic DNA of different genes at various concentrations were used as templates for PCR using reagents from TransGen and Company A. The amplification products were analyzed by 1.0% agarose gel electrophoresis.


## TransStart® FastPfu DNA Polymerase (AP221)

- Fast: 4 kb/min extension speed.
- High fidelity: 54 times higher fidelity than common Taq enzymes and 3 times higher than common Pfu enzymes.
- Long fragment amplification: Up to 15 kb of genomic DNA fragments and up to 20 kb of Plasmid DNA fragments can be amplified.
- High specificity: Using the "TransStart" double-blocking method novel hot-start technology, primers and templates are simultaneously blocked.
- Strong amplification ability: the unique PCR Stimulant enhances the enzyme's ability to amplify complex templates.
- The amplified product is blunt-ended and can be directly cloned into pEASY®-Blunt series vectors.

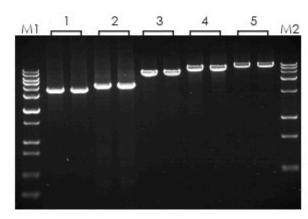

M1: 1Kb Plus DNA Ladder
M2: Trans15K DNA Marker
1: ACTR 3.5 kb;
2: VIN 4.6 kb;
3: Pol 6.8 kb;
4: APC 8.5 kb;
5: Dynein 12.3 kb;
cDNA: Amplification of cDNA with

TransStart FastPfu DNA Polymerase

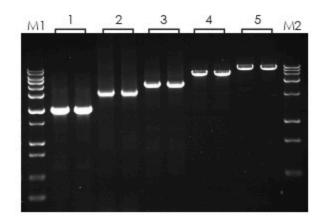

M1: 1Kb Plus DNA Ladder

M2: Trans 15K DNA Marker


1: Rhod 2.0 kb;

2: β-globin 3.0 kb;

3: Rhod 4.17 kb;

4: Factor IX 7.5 kb: 5: Serum albumin 12.4 kb;

Genomic DNA: Amplification of Genomic DNA with TransStart FastPfu DNA Polymerase

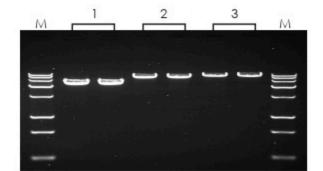

M: Trans15K DNA Marker

1: UDG 7.0 kb; 2: LN 10.0 kb;

3: Fang 14.7 kb;

Plasmid DNA: Amplification of Plasmid DNA with TransStart FastPfu DNA Polymerase


## TransStart® FastPfu Fly DNA Polymerase (AP231)

- Fast: 6 kb/min extension speed, 1kb extension takes only 10 seconds.
- High fidelity: 108 times higher fidelity than common Taq enzymes and 6 times higher than common Pfu enzymes.
- Long fragment amplification: Up to 15 kb of genomic DNA fragments and up to 20 kb of Plasmid DNA fragments can be amplified.
- High specificity: Using the "TransStart" double-blocking method novel hot-start technology, primers and templates are simultaneously blocked.
- Strong amplification ability: the unique PCR Stimulant enhances the enzyme's ability to amplify complex templates.
- The amplified product is blunt-ended and can be directly cloned into pEASY®-Blunt series vectors.

## **Bst II DNA Polymerase**


## (LP301)

- Optimized Buffer and enzyme SuperMix for Loop-Mediated Isothermal DNA Amplification (LAMP). Just add template and primers.
- Visualization: amplification results can be visualized by the color changing.
- Compatible with dUTP/UDG: good resistance to dUTP, adding dUTP/UDG to effectively prevent LAMP product contamination.

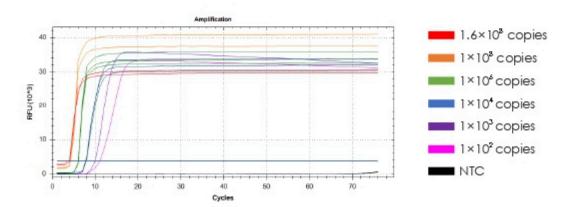

LAMP Amplification results from plasmid containing muscovy duck parvovirus (MDPV) as templats using TransGen Bst II DNA Polymerase.



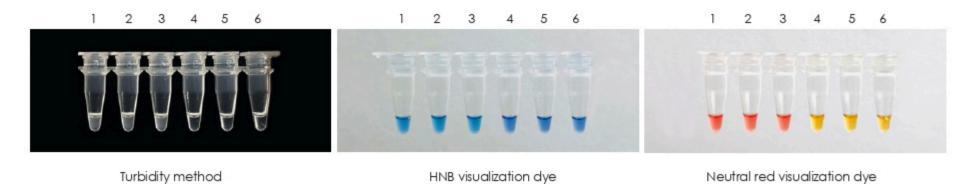


1-3: Plasmid containing ASFV-P72 gene 4-6: No Template Control (NTC)


## **Bst III DNA Polymerase** (LP311)

- Optimized Buffer and enzyme SuperMix for Loop-Mediated Isothermal DNA Amplification (LAMP). Just add template and primers.
- High efficiency: one-step method for reverse transcription and LAMP amplification.
- · Lyophilized format is available.



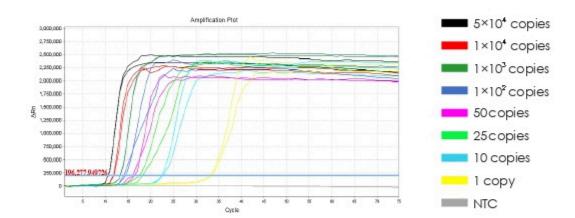


Different concentrations of in vitro transcribed duck and goose reovirus (DuRV) RNA were used as template for RT-LAMP amplification using TransGen Bst III DNA polymerase.

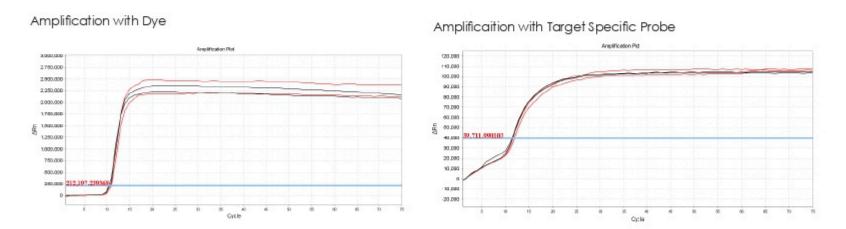

In vitro transcribed Duck and goose reovirus (DuRV) RNA was used as template for RT-LAMP amplification using TransGen Bst III DNA polymerase.

## Lyophilized Format vs Traditional Format

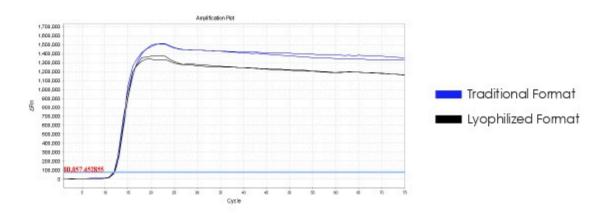


In vitro transcribed Duck and goose reovirus (DuRV) RNA was used as template for RT-LAMP amplification using conventional and lyophilized TransGen Bst III DNA polymerase. The results show that TransGen Bst III DNA polymerse can be lyophilized without lost its activity.


## T4 DNA Ligase for NGS (LL101)

- Optimized for adapter ligation during NGS library construction.
- Can be used for cloning of restriction enzyme digested nucleic acid fragments.



## Comparison with Competitive Product

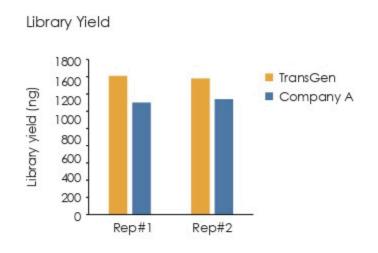

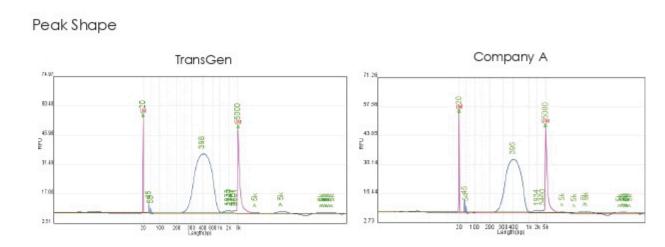

## DNA Polymerase I Klenow Fragment (LE201)

- Engineered from E.coli DNA polymerase that retains polymerase and 3'
   —> 5' exonuclease activity.
- Filling in at the 5' overhang of double-stranded DNA, and excision of the 3' overhang in double-stranded DNA to form blunt-ends.

## T4 DNA Polymerase (LP201)

- Catalyzes the synthesis of 3' to 5' DNA, and removes 3' overhangs or fills-in 5' overhangs to form blunt ends.
- Synthesis of labeled DNA probes by displacement reactions.

## T4 Polynucleotide Kinase (LK101)

- Catalyzes the phosphorylation of the 5' end of oligonucleotides, DNA or RNA.
- 5' end labeling of oligonucleotides, DNA or RNA. The labeled nucleic acids can used as probes for Southern blot, Northern blot, etc.





## Uracil-DNA Glycosylase (Low Temperature)

(LU201)


- Efficiency: the dU-containing template can be degraded as fast as 1 minute at room temperature.
- $\bullet$  Mininum inhibition to PCR reaction: No inhibition was observed when up to 4 Units of UDG was added to 50  $\mu$ l of PCR reaction.
- Wide range of applications: suitable for PCR/qPCR, RT-PCR/RT-qPCR, LAMP/RT-LAMP and other applications.



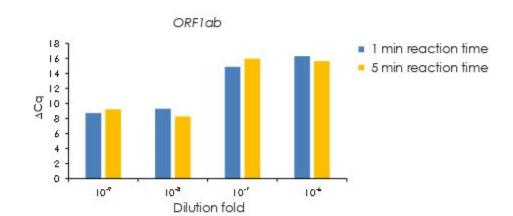

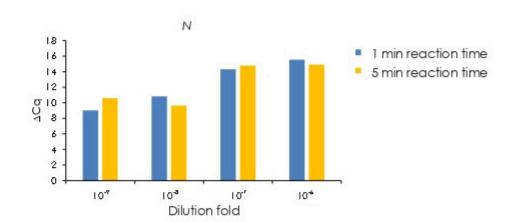


DNA fragments for ORFlab and N genes were amplified with PCR using dUTP to substitute dTTP. The amplified DNAs were diluted 10°, 10°, 10°, 10° folds, then treated with TransGen Uracil-DNA Glycosylase. The treated nucleic acids were used as templates for qPCR.

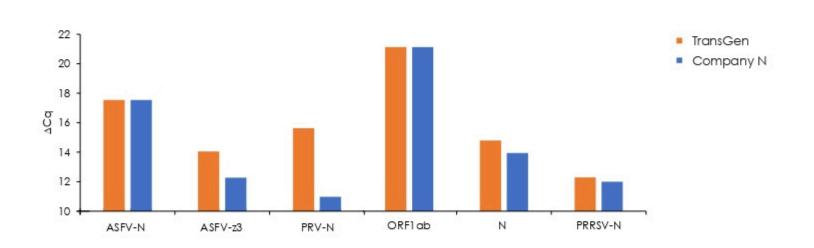

DNA fragments were amplified by PCR using dUTP to substitute dTTP. The amplified DNAs were treated with Uracil-DNA Glycosylase from TranGen or from Company N. The treated nucleic acids were used as templates to qualtify six genes by qPCR. (\DCQ represents the Cq difference of the three reaction systems of TransGen 2.5 U/mL UDG, Company N 2.5 U/mL UDG and 0 UDG)


## Low Inhibition of PCR Reactions

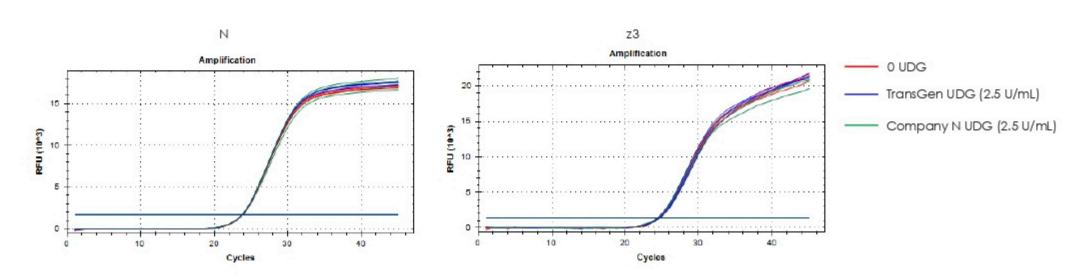




## Universal Nuclease (GMP Grade) (LN201)


- Purity ≥99%, specific activity ≥1.0×106 U/mg.
- Maintain the its activity under the following tested conditions: 6 M Urea, 0.1 M Guanidine HCl, 0.4% Triton X-100, 0.1% SDS, 1 mM EDTA, 1 mM PMSF.
- Compatible with various cell lysates (such as RIPA based) or protein extraction reagents.
- Compliant with pharmacopoeia requirements: no animal origin, no antibiotics, no endotoxin.
- Wide range of applications: degrades all forms of DNA and RNA and is widely used to remove nucleic acids from biological

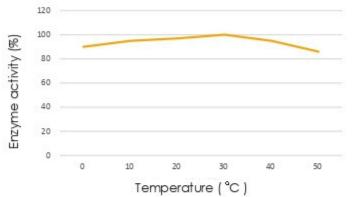


Universal Nuclease is active from temperature of 0-50°C.

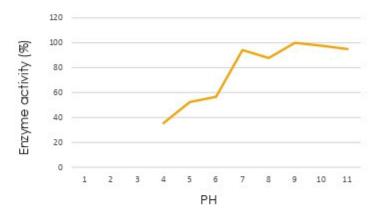

Universal Nuclease is active in the PH of 7-11.

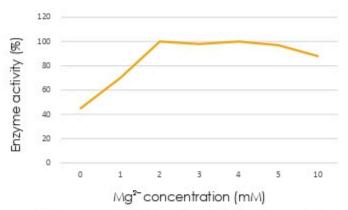

Universal Nuclease is active at Mg<sup>2-</sup> concentrations of 2-5 mM

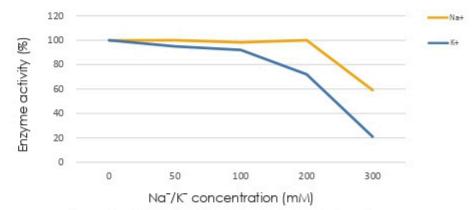

Universal Nuclease is active at Na¯ concentration of 0-200 mM and K¯ concentration of 0-100 mM.

## Quality Control

Item	Standard	Method		
Appearance	Colorless and transparent	Visual inspection		
Purity	≥99%	HPLC		
Enzyme activity	≥ 250U/µI	General substrate method		
Specific activity	≥1.0×106 U/mg	Enzymatic activity/protein concentration		
Bacteria residue	Not detectable	Culture method		
Protease activity	Not detectable	General substrate method		
Endotoxin content	< 0.25 EU/1000 U	Gel clo		
Host protein residue	≤10 ppm (µg/mL)	ELISA		
Mycoplasma residue	Not detectable	Mycoplasma detection kit (qPCR method)		
Pathogen detection (HCV/HBV/HIV)	Not detectable	PCR		

## TRANSGEN


## RT-PCR



## 01


## TransScript® IV One-Step gDNA Removal and cDNA Synthesis SuperMix (AW311)

- One reaction for gDNA remove and RT: Unique gDNA Remover for simultaneous reverse transcription and genomic DNA removal.
- 5 seconds to inactive gDNA remover and RT: gDNA Remover and reverse transcriptase are simultaneously heat-inactivated at 85°C for 5 seconds.
- High thermal stability: reverse transcription reaction temperature at 42°C - 65°C.
- High synthesis ability: up to 20 kb cDNA in length.

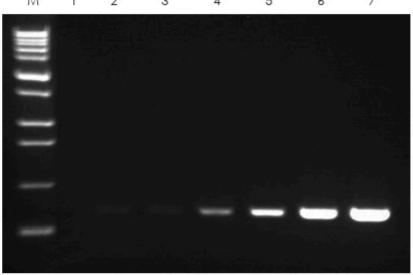

M: Trans2K<sup>®</sup> Plus II DNA Marker Total RNA used for RT-PCR: lane #1, 1 pg; lane #2,10 pg; lane #3, 100 pg; lane #4, 1ng; lane #5, 10 ng; lane #6,100 ng; lane #7, 1 µa.

RT-PCR were performed using *TransScript* IV OneStep gDNA Removal and cDNA Synthesis SuperMix with different concentrations of human total RNA as templates. The amplified nucleic acids were analyzed by agarose gel electrophoresis.

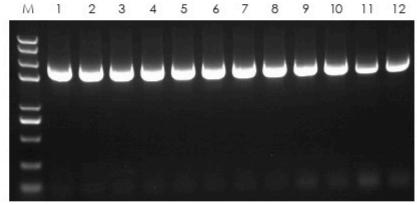

M: Trans2K® Plus II DNA Marker Amplified gene: BACH1 1-2: 50 °C 3-4: 55 °C 5-6: 56.4 °C 7-8: 60 °C 9-10: 62.5 °C 11-12: 65 °C


RT-PCR was performed with 100 ng of human total RNA at different temperature.

## TransScript® One-Step gDNA Removal and cDNA Synthesis SuperMix

(AT311)

- One reaction for gDNA remove and RT: Unique gDNA Remover for simultaneous reverse transcription and genomic DNA removal.
- 5 seconds to inactive gDNA remover and RT: gDNA Remover and reverse transcriptase are simultaneously heat-inactivated at 85°C for 5 seconds.
- Thermal stability: reverse transcriptase reaction temperature at 42°C.
- Up to 12 kb cDNA in length.


## **Comparsion with Competitive Products**

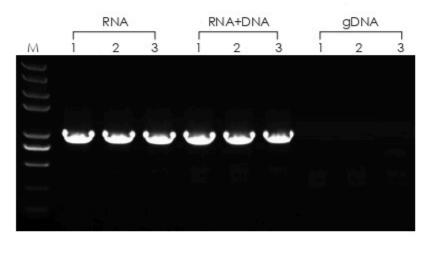

M: Trans2K® Plus II DNA Marker

40 35

- 1: TransGen 2: Company A
- 3: Company B

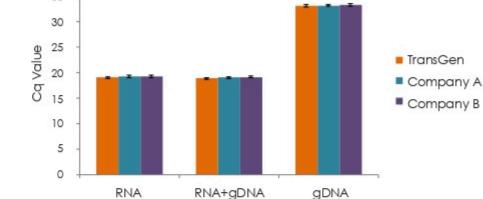


RT-PCR was performed with (1), 100 ng of human total RNA, (2), 100 ng of human total RNA and 200 ng gDNA, (3), 200 ng gDNA, using reagents from TransGen, Company A and Company B.

## TransScript® miRNA First-Strand cDNA Synthesis SuperMix (AT351)

- Optimized buffer for poly(A) polymerase and reverse transcriptase in optimized buffer, ensuring overall efficiency.
- Poly(A) tailing and cDNA synthesis in one reaction.
- Applicable for miRNA detection.

Comparsion with Competitive Product

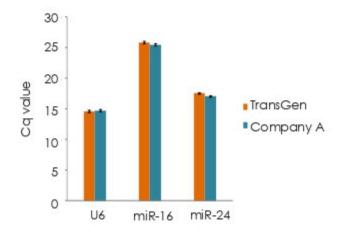

qRT-PCR was performed for miRNA extracted from human plasma using reagens from TransGen and Company A respectively.


## TRANSGEN qPCR&qRT-PCR





TransScript® II Multiplex Probe One-Step qRT-PCR SuperMix UDG

(AQ322)


- UDGase and dUTP to effectively prevent PCR contamination.
- Optimized for multiplex probe based qRT-PCR.
- High stability: reagents have been tested for repeated freezing and thawing for 30 cycles,
   stored at room temperature for 15 days, stored at 37°C for 7 days.
- Lyophilized format is also available.

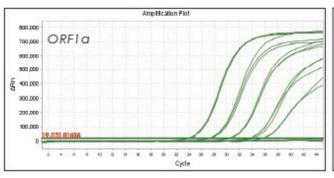


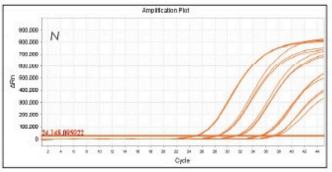


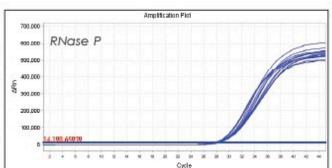


A serial of dilution of in vitro transcribed RNAs (5×10<sup>6</sup> copies/mL to 5×10<sup>2</sup> copies/mL, 10-fold dilution) were used as templates for qRT-PCR using TransScript II Multiplex Probe One-Step qRT-PCR Super/Mix UDG.

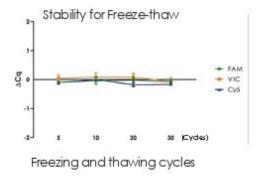

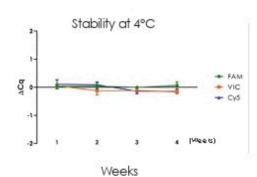

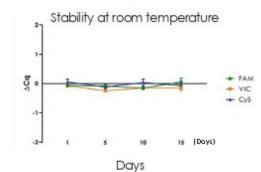

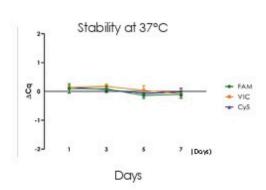



## PerfectStart® II Probe qPCR SuperMix UDG (AQ712)

- UDGase and dUTP to effectively prevent PCR contamination.
- Three different antibodies used for high specificity, high sensitivity, and high amplification efficiency.
- Optimized for probe based qPCR.
- High stability: reagents have been tested for repeated freezing and thawing for 20 cycles.
- Lyophilized format is also available.




## **Comparsion with Competitive Products**

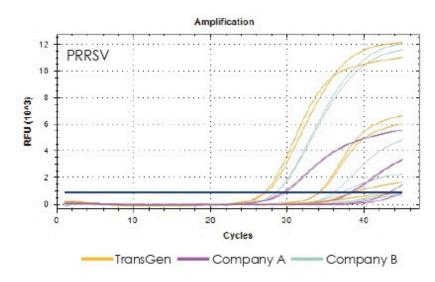

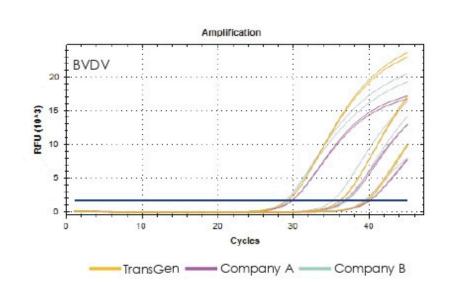

Different concentrations of cDNA (10 pg, 0.1 pg, 0.01 pg) from porcine reproductive and respiratory syndrome virus (PRRSV) and cDNA (100 pg, 1 pg, 0.1 pg) from bovine viral diarrhea virus (BVDV) were used as templates to perform probe based qPCR using reagents from TransGen, Company A, and Company B.

## Freeze-Thaw Study (Quadruplex)

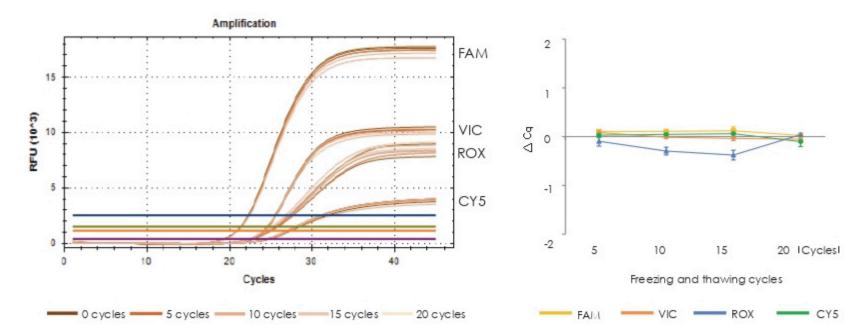


Quadruplex qPCR assays were performed using PerfectStart® II Probe qPCR SuperMix UDG after 5, 10, 15, and 20 cycles of freeze-thaw.


## PerfectStart® Fast Green qPCR SuperMix (AQ611)



- Three antibodies used for high specificity, high sensitivity, and high amplification efficiency.
- Fast amplification with 15 seconds extension time.
- Enhance specificity and reduce primer dimer formation with double cationic buffer.




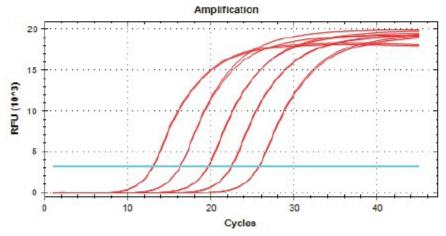


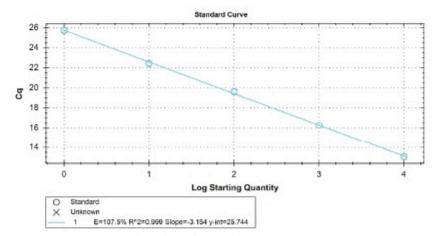

With 10 fold diluted plasmid DNA (1 ng to 0.1 pg) as template, qPCR was performed with PerfectStart® Fast Green qPCR SuperMix.


Standard two-step qPCR:

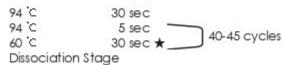

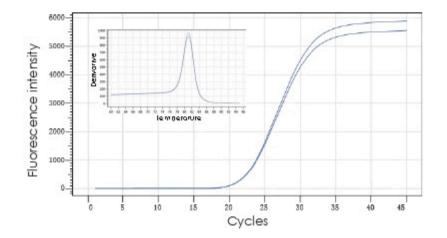


Two-step qPCR with PerfectStart® Fast Green qPCR SuperMix:



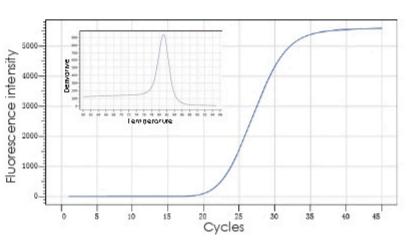

## TRANSGEN

## In Vitro mRNA Synthesis

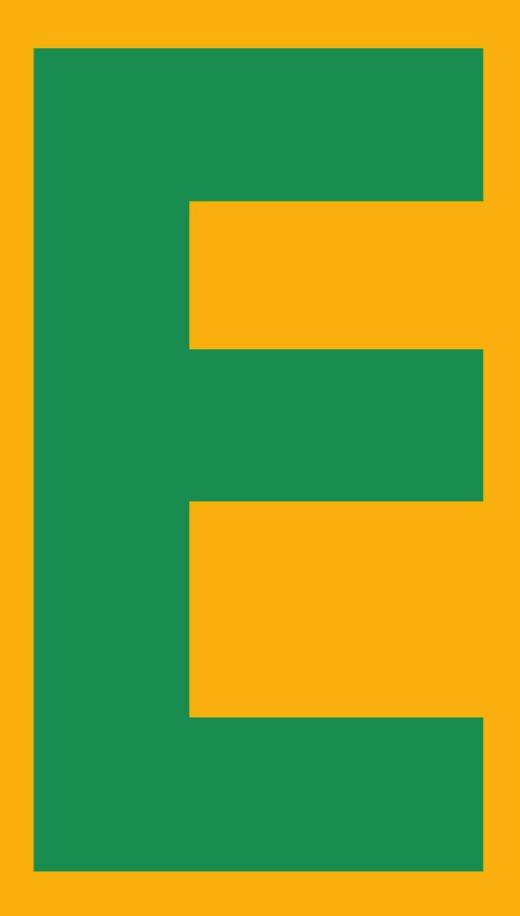

## TransGen In Vitro mRNA Synthesis Related Products



· EasyPure® HiPure PlasmidMaxiPrep/ MiniPrep Kit Phi29 DNA Polymerase ·FlyCut®Series Fast

**Endonucleases** 

- · T7 High Efficiency Transcription Kit · Ribonuclease Inhibitor
- ·DNase I (RNase-free) · Ppase


· mRNA Poly(A) Tailing


- · mRNA Capping Kit
- · MagicPure® Size **Selection DNA Beads** · MagicPure® RNA Beads
- · MagicPure®mRNA Kit
- · EasyPure® RNA **Purification Kit**
- · PerfectStart® E.coli DNA **Quantification Kit**

## **T7 High Efficiency Transcription Kit** (JT101)

- High yields: more than 200 µg of RNA can be produced in a 20 µL reaction volume, and milligrams of mRNA can be prepared in a single 1 mL reaction.
- Good for a variety of templates and different amount of templates: templates can be plasmids, PCR purified products, and synthetic DNA. Different amount of templates (1 ng to 2 µg) can be used for each reaction. Up to 6 kb of the transcription fragments can be obtained.

## RNA yield from 2 hours reaction time

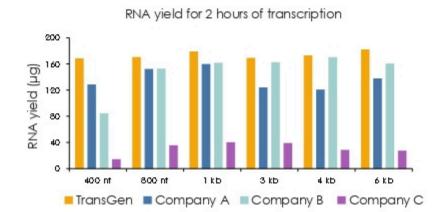


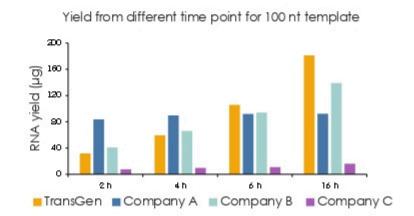


## 

## mRNA Capping Kit (LC101)

Adding vaccinia virus capping enzyme and cap structure 2'-O-methyltransferase to the target mRNA in the same reaction.
 Cap1 capped products can be obtained within 1 hour.

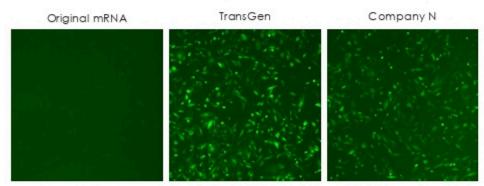

Fluorescent microscopy 48 hours post transfection

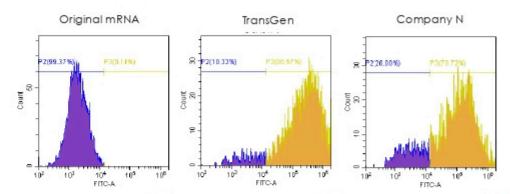


Transfection efficiency:0.14% Transfection efficiency:86.97% Transfection efficiency:79.72% Flow cytometric sorting 48 hours post transfection

## 

## mRNA Poly(A) Tailing Kit

- Adds a poly(A) tail of at least 150 nucleotides in length to the 3' termini of RNA.
- Reaction time can be adjusted to control for tail length.
- Capped mRNA from TransGen mRNA Capping Kit (LA201) can be directly used for the tailing reaction.
- Suitable for the polyadenylation of in vitro transcribed RNA to enhance translation efficiency.

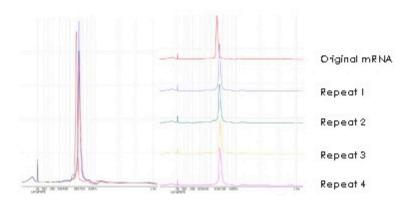

Classification	Product Name	Cat. No.	Specification
Nucleic Acid Purification	MagicPure® Blood Genomic DNA Kit (with Magnetic Stand)	EC101-01	50 rxns
	MagicPure® Blood Genomic DNA Kit (without Magnetic Stand)	EC101-11	50 rxns
	MagicPure® 32 Blood Genomic DNA Kit	EC101-32-11/12	32/64 rxns
	MagicPure® 96 Blood Genomic DNA Kit	EC101-96-11	96 rxns
	MagicPure® Fungi Genomic DNA Kit  MagicPure® 32 Fungi Genomic DNA Kit	EC103-01 EC103-32-11/12	50 rxns 32/64 rxns
	MagicPure® Cell-Free DNA Kit II (with Magnetic Stand)	EC103-32-11/12	32/64 rxris 50 rxns
	MagicPure® Cell-Free DNA Kit II (without Magnetic Stand)	EC211-01	50 rxns
	MagicPure® Viral DNA/RNA Kit (with Magnetic Stand)	EC301-01	50 rxns
	MagicPure® Viral DNA/RNA Kit (without Magnetic Stand)	EC301-11	50 rxns
	MagicPure® 32 Viral DNA/RNA Kit	EC301-32-11/12	32/64 rxns
	MagicPure® 96 Viral DNA/RNA Kit	EC301-96-11	96 rxns
	MagicPure® Simple Viral DNA/RNA Kit (with Magnetic Stand)	EC311-01	50 rxns
	MagicPure® Simple Viral DNA/RNA Kit (without Magnetic Stand)	EC311-11	50 rxns
	MagicPure® Simple 32 Viral DNA/RNA Kit	EC311-32-11/12	32/64 rxns
	MagicPure® Simple 96 Viral DNA/RNA Kit	EC311-96-11	96 rxns
	MagicPure® FFPE Tissue Genomic DNA Kit	EC701-01	50 rxns
	MagicPure® Stool and Soil Genomic DNA Kit	EC801-11	50 rxns
	EasyPure® Stool Genomic DNA Kit	EE301-01	50 rxns
	EasyPure® Micro Genomic DNA Kit	EE181-01	50 rxns
-	EasyPure® FFPE Tissue Genomic DNA Kit  EasyPure® Buccal Swab Genomic DNA Kit	EE191-01 EE201-01	50 rxns 50 rxns
-	EasyPure® Viral DNA/RNA Kit	ER201-01/02	50/200 rxns
-	MagicPure® Host Cell Residual DNA Kit (with Magnetic Stand)	EH201-01	50 rxns
-	MagicPure® Host Cell Residual DNA Kit (without Magnetic Stand)	EH201-11	50 rxns
-	EasyPure® miRNA Kit	ER601-01	50 rxns
<del></del>	TransStart® TopTag DNA Polymerase	AP151-01/02/03	250/500/6×500 units
İ	TransStart® TopTaq DNA Polymerase (with 2.5 mM dNTPs)	AP151-11/12/13	250/500/6×500 units
Ì	PerfectStart® Tag DNA Polymerase	AP401-01/02/03	250/500/6×500 units
	PerfectStart® Taq DNA Polymerase (with 2.5 mM dNTPs)	AP401-11/12/13	250/500/6×500 units
	TransStart® FastPfu DNA Polymerase	AP221-01/02/03	250/500/6×500 units
[	TransStart® FastPfu DNA Polymerase (with 2.5 mM dNTPs)	AP221-11/12/13	250/500/6×500 units
	TransStart® FastPfu Fly DNA Polymerase	AP231-21/22/23	250/500/6×500 units
Common Used	Phi29 DNA Polymerase	LP101-01/02	250/5×250 units
Enzymes	Bst II DNA Polymerase	LP301-01/02	100/200 rxns
	Bst III DNA Polymerase	LP311-01/02	100/200 rxns
	Universal Nuclease (GMP Grade)	LN201-01/02	25/50 KU
	Uracil-DNA Glycosylase (Low Temperature)	LU201-01/02	100/500 units
	T4 DNA Ligase(for NGS) T4 DNA Polymerase	LL101-01/02 LP201-01/02	200 µL / 1 mL 150 units / 5×150 units
	T4 Polynucleotide Kinase	LK101-01/02	500 units / 4×500 units
-	DNA Polymerase I Klenow Fragment	LE201-01	500 Units
	TransScript®    Reverse Transcriptase[M-MLV,RNaseH] (High Temperature RT)	AH101-02	10000 units
	TransScript® One-Step gDNA Removal and cDNA Synthesis SuperMix	AT311-02/03/04	50/100/500 rxns×20 µl reaction
RT-PCR -	TransScript® One-Step RT-PCR SuperMix	AT411-02	200 rxns×20 µl reaction
	TransScript® Uni One-Step RT-PCR SuperMix	AU411-02	200 rxns×20 µl reaction
	TransScript® IV One-Step gDNA Removal and cDNA Synthesis SuperMix	AW311-02/03	50/100 rxns
	TransScript® miRNA First-Strand cDNA Synthesis SuperMix	AT351-01	20 rxns×20 µl reaction
	TransScript® Green miRNA Two-Step qRT-PCR SuperMix	AQ202-01	RT reaction/qPCR reaction 20 rxns×20 µl/500 rxns×20 µl
	TransScript® Green One-Step qRT-PCR SuperMix	AQ211-01/02	100/400 rxns×20 µl reaction
	PerfectStart® II Probe qPCR SuperMix UDG	AQ712-01/02/03	1/5×1/15×1 ml
qPCR & qRT-PCR	TransScript® II Multiplex Probe One-Step qRT-PCR SuperMix UDG	AQ322-01/02	100/400 rxns×20 µl reaction
	PerfectStart® Green qPCR SuperMix	AQ601-01/02/03/04	1/5×1/15×1/25×1 ml
	PerfectStart® Fast Green qPCR SuperMix	AQ611-01/02/03/04	1/5×1/15×1/25×1 ml
	TransScript® Probe One-Step qRT-PCR SuperMix	AQ221-01/02	100/400 rxns×20 µl reaction
In Vitro mRNA Synthesis	EasyPure® HiPure Plasmid MiniPrep Kit	EM111-01	50 rxns
	EasyPure® HiPure Plasmid MaxiPrep Kit	EM121-01-V2	10 rxns
	FlyCut® Fast Endonuclease	See catalog	See catalog
	T7 High Efficiency Transcription Kit	JT101-01/02	25/100 rxns×20 μl
	mRNA Capping Kit	LC101-01/02	25/100 rxns
	mRNA Poly(A) Tailing Kit	LA201-01/02	25/100 rxns
	MagicPure® Size Selection DNA Beads	EC401-01/02/03/04	1/5/60/450 ml
	MagicPure® RNA Beads MagicPure® mRNA Kit	EC501-01/02/03	1/5/60 ml
	MagicPure® RNA Purification Kit	EC511-01/02 ER701-01	24/96 rxns 25 rxns
	RNAhold®	EH101-01	25 rxns 100 ml
ŀ	High Pure dNTPs (2.5 mM)	AD101-01/02	1/5×1 ml
Delete d	High Pure dNTPs (10 mM)	AD101-01/02 AD101-11/12	1/5×1 ml
Related	and the second control of the second control	AI101-01/02	2000/5×2000 units
76 CCC CC	Riponuciease innibitor		ern arriver and a comment of the com
Products	Ribonuclease Inhibitor  DNase I (RNase-free)		
76 CCC CC	DNase I (RNase-free)  Proteinase K	GD201-01 GE201-01	1 500 units

## CliniSciences Group

## Austria

Company: CliniSciences GmbH Address: Sternwartestrasse 76, A-1180

Wien - Austria

Telephone: +43 720 115 580 Fax: +43 720 115 577

Email: <u>oesterreich@clinisciences.com</u>
Web: <u>https://www.clinisciences.com</u>



Belgium Company: CliniSciences S.R.L Address: Avenue Stalingrad 52, 1000

Brussels - Belgium

Telephone: +32 2 31 50 800 Fax: +32 2 31 50 801

Email: belgium@clinisciences.com Web: https://www.clinisciences.com



## Denmark

Company: CliniSciences ApS Address: Oesterbrogade 226, st. 1, Copenhagen, 2100 - Denmark Telephone: +45 89 888 349 Fax: +45 89 884 064

Email: danmark@clinisciences.com Web: https://www.clinisciences.com



## Finland

Company: CliniSciences ApS Address: Oesterbrogade 226, st. 1, Copenhagen, 2100 - Denmark Telephone: +45 89 888 349 Fax: +45 89 884 064

Email: suomi@clinisciences.com Web: https://www.clinisciences.com



## France

Company: CliniSciences S.A.S Address: 74 Rue des Suisses, 92000

Nanterre- France

Telephone: +33 9 77 40 09 09
Fax: +33 9 77 40 10 11
Email: info@clinisciences.com
Web: https://www.clinisciences.com



## Germany

Company: Biotrend Chemikalien GmbH Address: Wilhelm-Mauser-Str. 41-43,

50827 Köln - Germany
Telephone: +49 221 9498 320
Fax: +49 221 9498 325
Email: info@biotrend.com
Web: https://www.biotrend.com



## Iceland

Company: CliniSciences ApS Address: Oesterbrogade 226, st. 1, Copenhagen, 2100 - Denmark Telephone: +45 89 888 349 Fax: +45 89 884 064

Email: island@clinisciences.com Web: https://www.clinisciences.com



## Ireland

Norway

Company: CliniSciences Limited
Address: Ground Floor, 71 lower Baggot street

Dublin D02 P593 - Ireland Telephone: +353 1 6971 146 Fax: +353 1 6971 147

Company: CliniSciences ApS

Copenhagen, 2100 - Denmark

Telephone: +45 89 888 349

Fax: +45 89 884 064

Address: Oesterbrogade 226, st. 1,

Email: norge@clinisciences.com

Web: https://www.clinisciences.com

Email: <u>ireland@clinisciences.com</u>
Web: <u>https://www.clinisciences.com</u>



## Italy

Company: CliniSciences S.r.l Address: Via Maremmana inferiore 378 Roma 00012 Guidonia Montecelio - Italy

Web: https://www.clinisciences.com

Telephone: +39 06 94 80 56 71 Fax: +39 06 94 80 00 21 Email: italia@clinisciences.com



## Netherlands

Company: CliniSciences B.V. Address: Kraijenhoffstraat 137A, 1018RG Amsterdam, Netherlands Telephone: +31 85 2082 351 Fax: +31 85 2082 353

Email: nederland@clinisciences.com Web: https://www.clinisciences.com



## Dolond

Company: CliniSciences sp.Z.o.o.
Address: ul. Rotmistrza Witolda Pileckiego 67
lok. 200 - 02-781 Warszawa -Poland

Telephone: +48 22 307 0535 Fax: +48 22 307 0532

Email: polska@clinisciences.com Web: https://www.clinisciences.com



## Portugal

Company: Quimigen Unipessoal LDA Address: Rua Almada Negreiros, Lote 5, Loja 14, 2615-275 Alverca Do Ribatejo - Portugal

Telephone: +351 30 8808 050
Fax: +351 30 8808 052
Email: info@quimigen.com
Web: https://www.quimigen.pt



Spain
Company: CliniSciences Lab Solutions
Address: C/ Hermanos del Moral 13

Address: C/ Hermanos del Moral 13 (Bajo E), 28019, Madrid - Spain Telephone: +34 91 269 40 65 Fax: +34 91 269 40 74

Email: <a href="mailto:espana@clinisciences.com">espana@clinisciences.com</a>
Web: <a href="mailto:https://www.clinisciences.com">https://www.clinisciences.com</a>



## Sweden

Company: CliniSciences ApS Address: Oesterbrogade 226, st. 1, Copenhagen, 2100 - Denmark Telephone: +45 89 888 349 Fax: +45 89 884 064

Email: sverige@clinisciences.com
Web: https://www.clinisciences.com



## Switzerland

Company: CliniSciences Limited Address: Marktgasse 18 8302 Kloten -

Switzerland

Telephone: +41 (044) 805 76 81 Fax: +41 (044) 805 76 75

Email: switzerland@clinisciences.com
Web: https://www.clinisciences.com



## UK

Company: CliniSciences Limited
Address: 11 Progress Business center, Whittle
Parkway, SL1 6DQ Slough- United Kingdom
Telephone: +44 (0)1753 866 511

or +44 (0) 330 684 0982 Fax: +44 (0)1753 208 899 Email: <u>uk@clinisciences.com</u> Web: <u>https://www.clinisciences.com</u>



Company: Biotrend Chemicals LLC Address: c/o Carr Riggs Ingram, 500 Grand Boulevard, Suite 210 Miramar

Beach, FL 32550- USA
Telephone: +1 850 650 7790
Fax: +1 850 650 4383
Email: info@biotrend-usa.com
Web: https://www.biotrend-usa.com

